

Endoluminal stentgraft system for open heart surgery

Ordering information


We deliver E-vita open hybrid stentgraft in diameters of 24-40 mm. Product training as well as on-site implantation support is part of our services.

Quick and reliable delivery of the standard products in Europe within 24 hours.

Call our central office: +49 (0) 74 71 / 922 - 0

E-vita open

Cuff length (mm)	Cuff ϕ (mm)	ϕ proximal Stentgraft (mm)	ϕ distal Stentgraft (mm)	Length Stentgraft (mm)	Catalogue No.
70	24	24	24	150	71AC2424S15-00
70	28	28	28	150	71AC2828S15-00
70	30	30	30	150	71AC3030S15-00
70	32	33	33	160	71AC3333S16-00
70	34	36	36	160	71AC3636S16-00
70	38	40	40	160	71AC4040S16-00

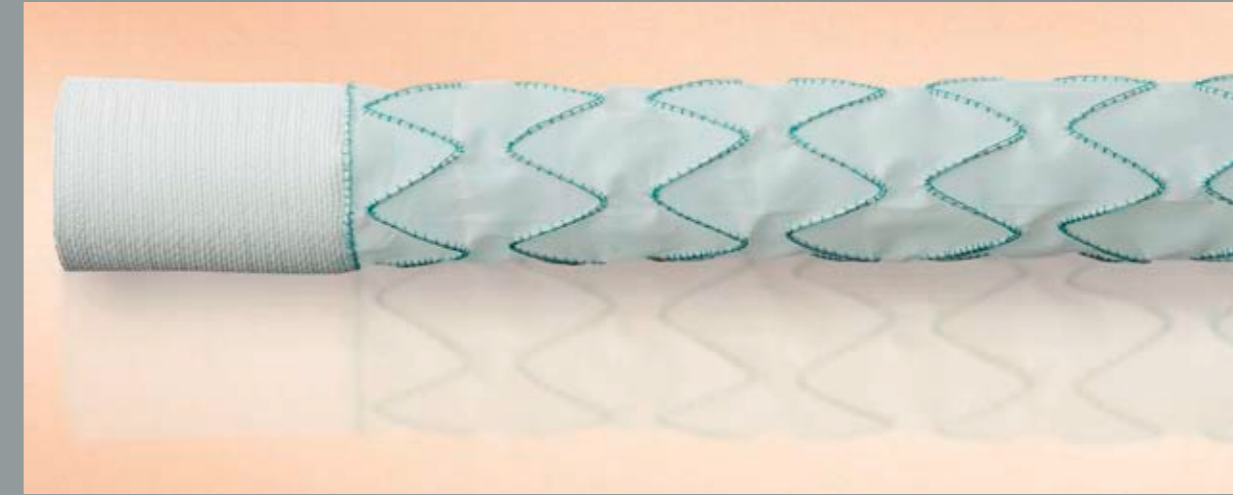

The endoluminal stentgraft system for open heart surgery

The hybrid system with clear advantages for patient and cardiac surgeon



Short and easy E-vita open delivery system with atraumatic ballpoint tip


Optimized “elephant trunk” procedure through combination of conventional and interventional techniques

The following indications can be treated in a one-step procedure with E-vita open developed by JOTEC:

- acute aortic dissection Stanford type A
- complex dissection of the aorta Stanford type B
- long-range aortic aneurysm
- chronic long-range dissection of the aorta

Two JOTEC products, the FlowWeave woven polyester graft and the E-vita thoracic stentgraft both with extensive clinical history have been combined to make this unique product: **E-vita open**

The established component combination of polyester for the prosthesis and Nitinol for the stent structure guarantees a high degree of product quality and safety.

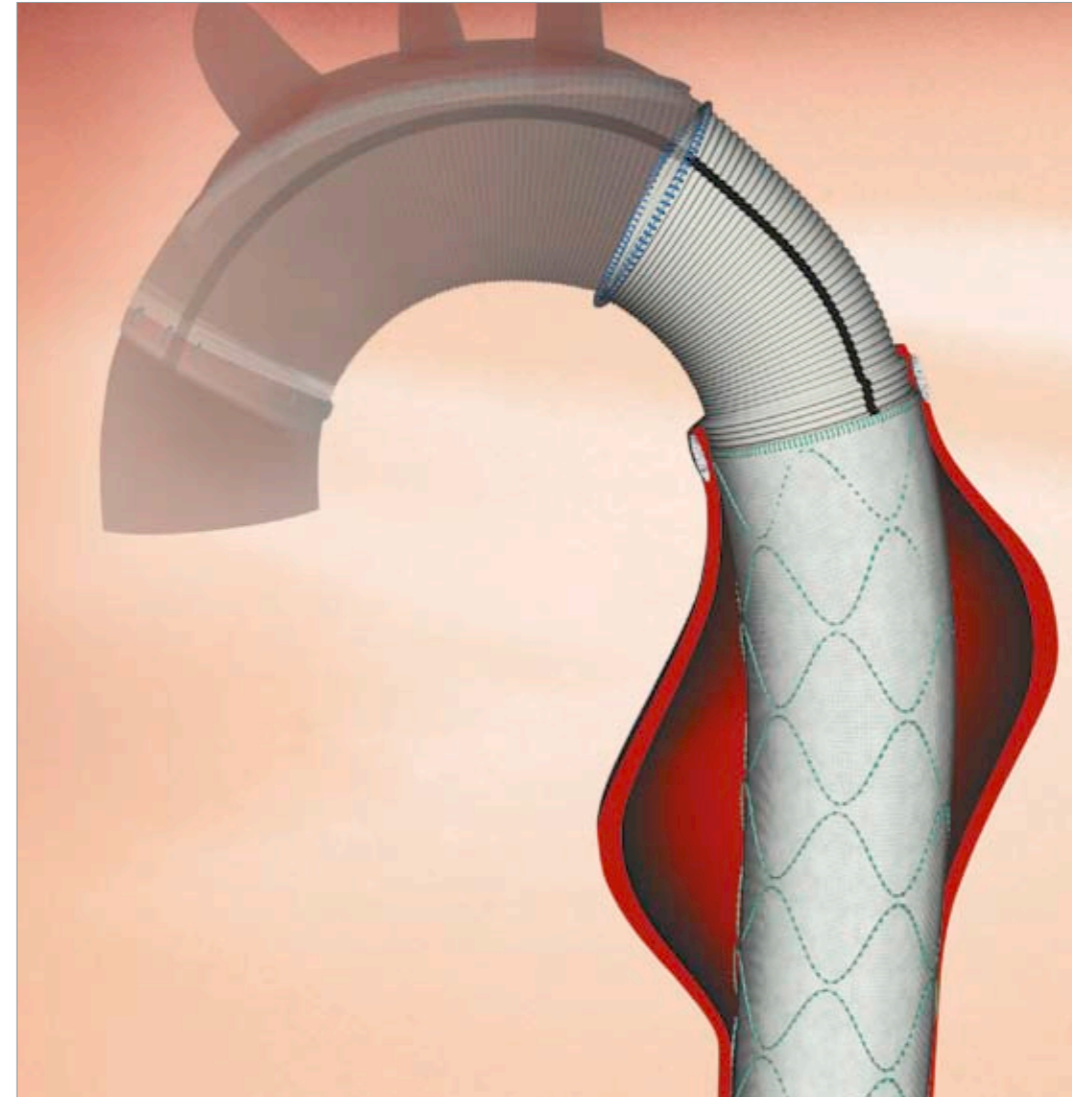

E-vita open is a stentgraft system for the cardiac surgeon

Prompt availability, ease of use, time saving as well as cost efficiency are clear advantages for the surgeon. The patient benefits from a safe treatment avoiding a second intervention.



Ultimate precision as well as severest quality control result in highest standards of quality



Proximal part of E-vita open with polyester cuff



Atraumatic ballpoint tip prevents injuries of the vessel wall



Outstanding flexibility inside the vessel

E-vita open is the only system where a stentgraft with a cuff is loaded in a delivery system. The system allows a one-step intervention by combining classic vascular reconstruction with modern aortic stenting significantly simplifying previous therapeutic procedures.

The stentgraft of the E-vita open is ideal when treating the non-accessible descending areas of the thoracic aorta. The woven prosthetic component allows a safe surgical fixation and serves as a connection (vascular adaptor) for the classical vascular reconstruction of the aortic arch.